ceramill therm s

DE Betriebsanleitung	3 - 44
EN User Manual	45 - 86
FR Mode d'emploi	87 - 130
IT Istruzioni d'uso	131 - 174
ES Modo de empleo	175 - 217

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Tab. 9

- Translation of the original Instruction Manual -

Table of Contents

1	Explanation of Symbols	46	10	Operation	78
2	General Safety Instructions	46	10.1	Switching on and off	78
3	Suitable Personnel	47	10.2	Placing on and removing a firing object	78
4	Machine Specifications	48	10.3	Selecting and starting programmes	79
4.1	Delivery Scope	48	10.4	Programme end	79
4.2	Intended Use	48	10.5	After a power failure	79
4.3	CE Declaration of Conformity	48	10.6	Access rights	80
4.4	Components and Interfaces	49	11	Cleaning and Maintenance	80
5	Installation	51	11.1	Cleaning	80
5.1	Setting up	51	11.2	Maintenance	81
5.2	Electrical connection	51	12	Storage	84
5.3	Initial starting of operation	52	13	Malfunctions, Repairs and Warranty ..	84
5.4	Connecting the furnace with the network (optional)	52	13.1	Malfunctions	84
5.5	Connecting the furnace with an individual computer (optional)	56	13.2	Repairs	84
6	Operating elements	58	13.3	Warranty	84
7	Screen displays	60	14	Environmental Protection	84
7.1	Programme Manager	61	15	Technical Data and Accessories/Spare Parts	85
7.2	Programme indication	64			
7.3	Programme sequence	66			
8	Programme editor	68			
8.1	CVTL commands	68			
8.2	Calling up the programme editor ..	69			
8.3	Cancelling the programme editor ..	69			
8.4	Creating a new programme	69			
9	Setup menu	70			
9.1	Settings	72			
9.2	Programmes	74			
9.3	Communication	75			

EN

1 Explanation of Symbols

Warning indications

Warning indications in the text are marked with a triangle and boxed.

In case of hazards through electricity, the exclamation mark in the warning triangle is substituted by a lightning bolt.

When there is a risk of burns from hot surfaces, the exclamation mark in the warning triangle is substituted by a heat-alert symbol.

Signal words at the beginning of a warning indication specify the type and severity of the consequences, if the measures to avert the hazard are not adhered to.

- _ **NOTE** means that property damage can occur.
- _ **CAUTION** means that light to fairly serious personal injury can occur.
- _ **WARNING** means that serious personal injury can occur.
- _ **DANGER** means that life-threatening personal injury can occur.

Important information

Important information that do not lead to hazards for humans or property damage are marked with the icon aside and are also boxed.

Other symbols in the Manual

Symbol	Meaning
▷	Item of an operation description
–	Item of a list
▪	Subitem of an operation description or a list
[3]	Numbers in square brackets refer to position numbers in graphics/figures

2 General Safety Instructions

When installing, starting-up and operating the machine, always observe the following safety instructions:

DANGER:

Hot surfaces

During operation, the firing chamber can reach temperatures above 1500 °C. There is acute danger of burning when touching hot surfaces after opening the firing chamber. After switching off, the furnace requires several hours to cool down.

- ▷ Under all circumstances, avoid direct contact with hot components of the furnace.
- ▷ Do not reach into the firing chamber with your hands or other objects.
- ▷ Do not touch the lift plate, insulation table or firing objects with your bare hands.
- ▷ Use suitable aids (tweezers) when removing firing objects or wear safety gloves.

**DANGER:**

This product does not comply with the ATEX Directive.

- ▷ Do not use in explosive atmospheres!
- ▷ Do not operate with explosive gases or mixtures!

**WARNING:**

Operation with power sources, products, consumables, auxiliary aids, etc., which are subject to the Ordinance on Hazardous Substances or in any way have an impact on the health of the operating personnel, is not permitted.

**WARNING:**

Fire hazard!

- ▷ When setting up the furnace, provide for suitable fire-protection measures (e. g., fire-resistant table and floor covering, fire-resistant protective walls).
- ▷ Never start the operation of the furnace without the insulation table mounted.
- ▷ Do not cover off the ventilation openings in the firing-chamber lid. Pay attention that combustible materials (e. g., rubbing alcohol, wipes, etc.) and heat-sensitive objects have sufficient clearance to the furnace.
- ▷ In case of unexpected processes in the furnace (for example, strong smoke development or odour nuisance), switch the furnace off immediately and allow the furnace to cool down to room temperature.

**NOTE:**

- ▷ Switch the machine off when not in use or unsupervised for longer periods, e.g., overnight. This method of saving electrical energy also benefits the environment.

3 Suitable Personnel

**NOTE:**

Starting-up and operation of the machine may only be carried out by trained specialised personnel.

**CAUTION:**

A defective machine can lead to malfunctions! On detection of damage or a functional defect of the machine:

- ▷ Label the machine as defective.
- ▷ Prevent further operation until the machine has been repaired.

4 Machine Specifications

4.1 Delivery Scope

- _ High-temperature furnace Ceramill Therm s
- _ Insulation table
- _ Sinter ring
- _ Sinter plate
- _ Furnace tweezers
- _ Calibration set for temperature calibration

▷ After unpacking, check the machine for completeness and possible transport damages. Please claim any transport damages immediately with your supplier.

4.2 Intended Use

Ceramill Therm s is a high-temperature furnace for the final sintering of dental oxide ceramics, particularly for Ceramill zirconium oxide blanks. The sintering furnace is pre-programmed by the manufacturer with one speed-sintering programme and six optimised standard programmes, which are matched to the Ceramill zirconium oxide frame material.

Any other or further-going use, such as the processing of other products than those intended for, as well as the handling of hazardous substances or materials hazardous to one's health, does not apply as under the intended use.

As to the materials being used in the furnace, it must be known whether they are aggressive to or even could lead to destruction of the insulation or of the heating rods.

It is prohibited to remove safety devices, bypass them or put them out of operation.

Unauthorised modifications/alterations of the machine shall also void the warranty.

The set-up instructions and the safety regulations must be followed, as otherwise the furnace is considered as not being used as intended for, and any and all claims against Amann Girrbach shall become void.

4.3 CE Declaration of Conformity

This product was designed and manufactured based on careful selection of the harmonised standards to be observed, as well as additional technical specifications. It thus corresponds with the state-of-the-art and ensures maximum safety. In terms of design and performance, this product complies with the European Directives and the supplementary national requirements. Conformity has been confirmed with the CE marking.

EC Guidelines

2006/95/EC	2004/108/EC
------------	-------------

Applied harmonised standards

DIN EN ISO 12100	DIN EN ISO 13732-1
DIN EN 61010-1	DIN EN 61000-6-1, DIN EN 61000-6-3

Tab. 1

The product's Declaration of Conformity can also be viewed on the manufacturer's website.

4.4 Components and Interfaces

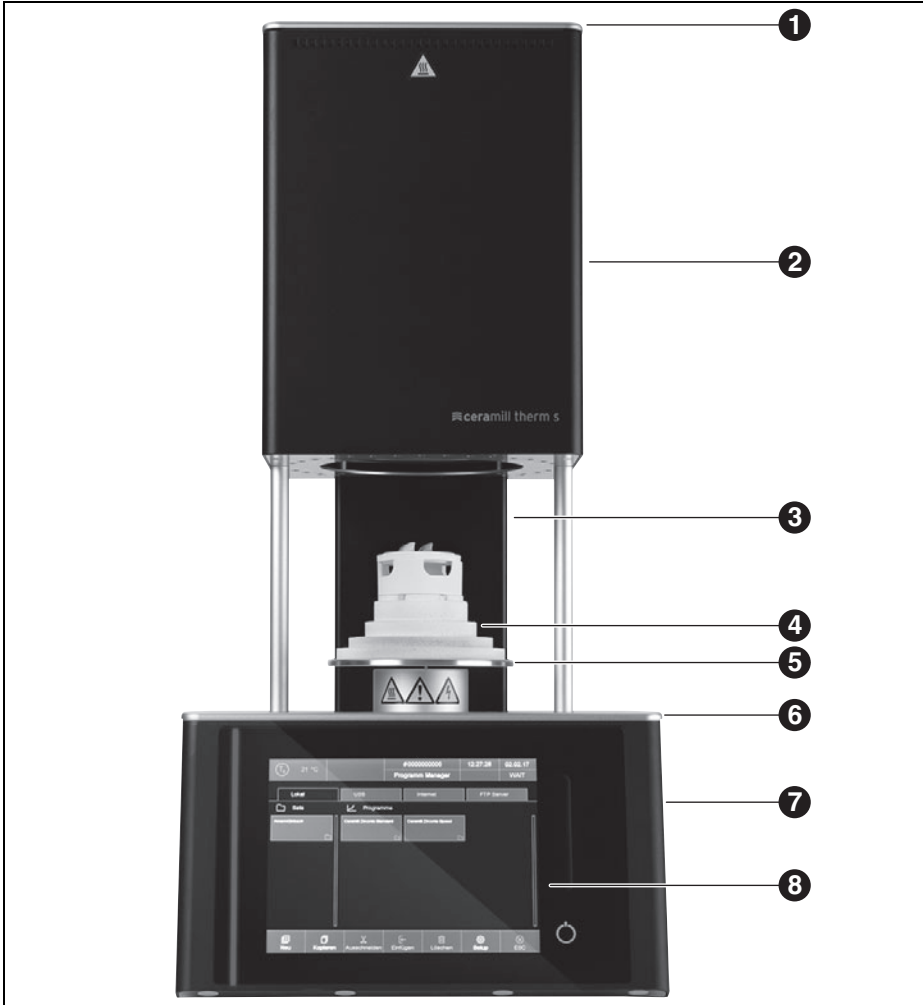
Fig. 1 Machine overview, front

- 1 Firing-chamber lid
- 2 Furnace top with firing chamber
- 3 Rear column
- 4 Insulation table
- 5 Lift plate with seal ring
- 6 Bottom cover plate
- 7 Furnace bottom with electronics
- 8 User interface with touchscreen

EN

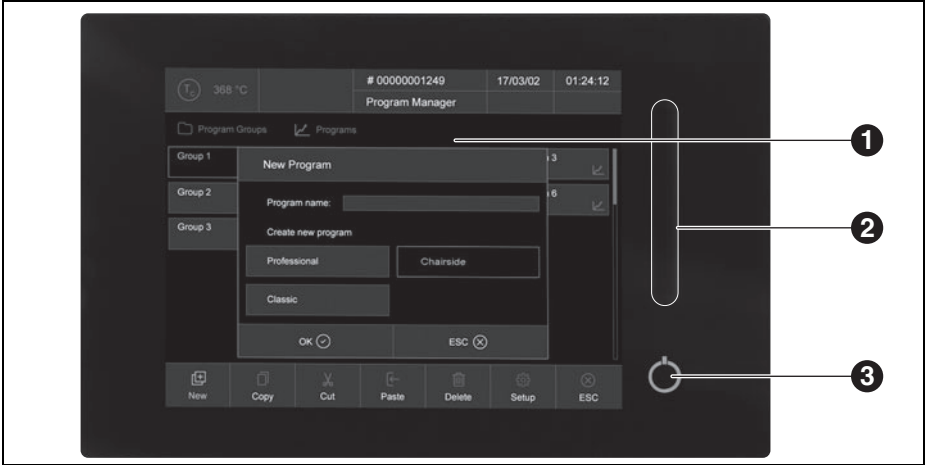
Fig. 2 User interface

- 1 Touchscreen
- 2 Swipe pad, lift
- 3 On/Off switch

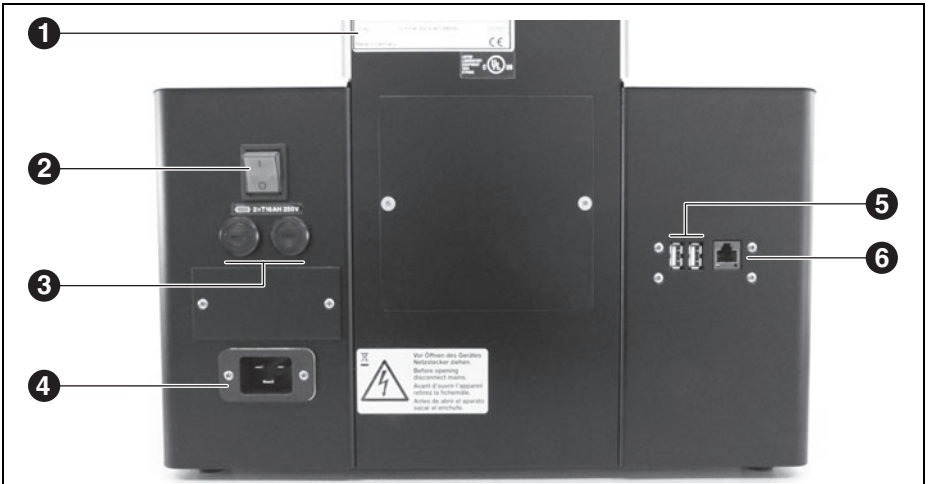
Fig. 3 Machine overview, rear

- 1 Type plate with serial number
- 2 Mains switch
- 3 Fuses
- 4 Power supply connection
- 5 USB ports
- 6 Ethernet connection

5 Installation

5.1 Setting up

 **NOTE:**
 Due to the weight of 23 kg, it is recommended that two persons set up the unit.
 ▷ When setting up, observe the valid country-specific accident prevention regulations.

- _ The machine is intended exclusively for use within dry, closed rooms.
- _ The machine must be positioned below an exhaust pipe or an exhaust hood (as used for pre-heating furnaces).
- _ A mains socket outlet must be located close to the unit and be easily accessible.
- _ Sufficient ventilation must be ensured at the set-up location, in order to draw off the heat and possibly developing exhaust gas.
- _ A clearance of at least 50 mm is to be observed sideways, to the rear and upwards of the furnace. The lateral clearance to non-flammable materials can be reduced to 20 cm. Space requirements for the Ceramill Therm s incl. connections is thus:
 - Width: 1.38 m (0.78 m for non-flammable materials)
 - Depth: 0.84 m
 - Height: 1.28 m
- _ The furnace ventilation opening must remain unobstructed to enable a cooling air circulation.
- _ The set-up surface must be level.
- _ The empty weight of the machine is 23 kg. The set-up surface must withstand an appropriate load.
- _ During operation, the room temperature must be between 15 °C and 35 °C; avoid large variations in temperature.

- _ The maximal operating elevation (standard elevation zero) is 2000 m.
- _ The floor covering must be made of non-flammable material, so that hot material falling out of the furnace can not inflame the covering.
- ▷ Set up the machine on a sturdy work table or work bench (on a non-flammable surface).
- ▷ Remove all packaging material, including from the furnace compartment.

 **DANGER:**
 For carrying out the further installation steps, a qualified professional in the area of electrical safety is mandatory. The installation and repair of heating elements must be performed by trained and certified personnel.

- ▷ Have the heating elements mounted by a trained professional (see separate installation instructions).

5.2 Electrical connection

The electrical connection is to be made via a socket outlet with earthing contact providing 220-240 V and 50 - 60 Hz. This socket outlet must be located close to the furnace and be easily accessible.

When using an extension cord or a multiple socket outlet, do not exceed the respective max. electrical load rating.

WARNING:

The mains cable may not be damaged!

- ▷ Do not place objects on the mains cable.
- ▷ Lay the mains cable in such a manner that nobody can step on or trip over it.
- ▷ Pay attention that the cable does not have contact with hot objects or surfaces.
- ▷ In case of damage to the insulation: Replace the mains cable.

NOTE:

Do not insert the mains plug before the heating elements have been inserted.

- ▷ Set the mains switch to the **0** position.
- ▷ When the furnace is to be operated via a network (Ethernet) or via a single computer: Connect network cable to the interface on the rear of the unit ([6] in Fig. 3 on page 50).

For a direct connection, older computers may require a crossover network cable.

- ▷ Connect the enclosed mains cable to the socket intended for this purpose on the rear of the furnace.
- ▷ Connect the other end of the cable to the socket outlet with earthing contact.

The mains cable may only be replaced against an original mains cable from the manufacturer.

5.3 Initial starting of operation

- ▷ Switch the furnace on via the mains switch.
The furnace control system boots up (approx. one minute). Afterwards, an internal check of the hardware and software takes place and the lift automatically moves to the bottom position.
- ▷ Check, if the insulation table is seated on the lift plate. If not: Switch the furnace off and mount the insulation table.

CAUTION:

Health hazard!

- ▷ Do not inhale ceramic dust.

- ▷ Clean ceramic dust from the lift plate and seal ring. Do not blow away the dust; carefully brush it into a receptacle using a brush.

The furnace's control system software is operated via the touchscreen.

5.4 Connecting the furnace with the network (optional)

5.4.1 Determining the IP address

If required, contact your network administrator.

- ▷ Make sure that the furnace is connected with the Ethernet.
- ▷ Switch the furnace on and confirm possible error messages.
After switching on, the Programme Manager is automatically displayed.

- ▷ In the Programme Manager, tap on *Setup > Netzwerk (Setup > Network)*.
The *Netzwerkeinstellungen (Network settings)* window opens.

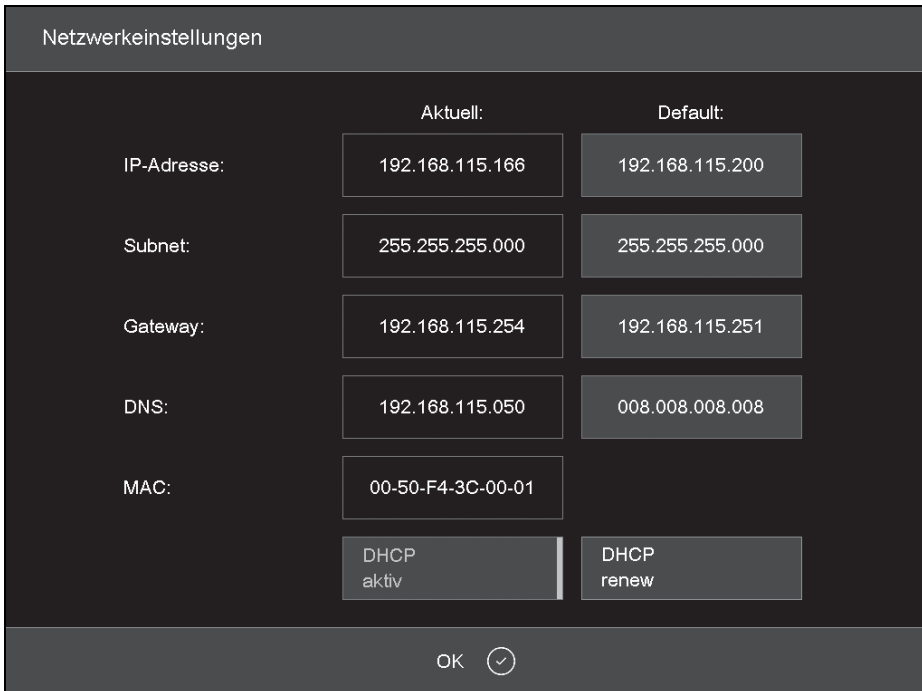
Fig. 4 *Setup > Network*

- ▷ Tap on the *DHCP aktiv (DHCP active)* button.
When the DHCP function is active, the button is highlighted in colour. When the DHCP function is active, the network server assigns the furnace an IP address.
- ▷ Tap on *OK*.
The message and the *Netzwerkeinstellungen (Network settings)* window are closed.
- ▷ Switch the furnace via the mains switch off and then on again.
This causes the control system to take over the changed network settings, and the furnace automatically seeks an IP address in the Ethernet network.
- ▷ In the Programme Manager, tap on *Setup > Netzwerk*.
The *Netzwerkeinstellungen (Network settings)* window opens.
- ▷ Note down the current IP address.

5.4.2 Setting up the furnace operation

Operation via browser

For operation of the furnace via Web interface, your browser must be Java-capable and Java must be activated.

- ▷ Start a Java-capable browser (e. g. Internet Explorer or Mozilla Firefox) on a computer in the network.
- ▷ Enter the noted down IP address in the browser's address line (e. g., "http://192.168.115.156") and press *Enter*. The Login window is displayed.

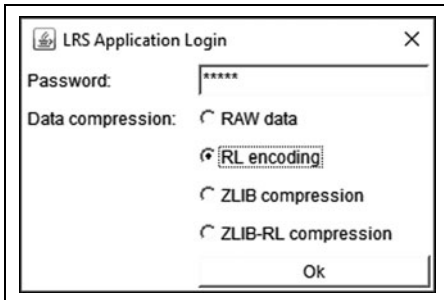
Fig. 5 Login window

- ▷ Enter the password in the *Password* entry line.
- ▷ Select the desired data compression for the data transfer between the furnace and the computer.
We recommend the standard setting *RL encoding*.
- ▷ Click *OK*.
The Programme Manager is displayed in the browser.
- ▷ Installing a fixed IP address for the furnace in the Programme Manager (see Chapter 5.4.3.).

For comfortable calling up of the furnace control system, you can define a bookmark in the browser:

- ▷ Enter the fixed IP address in the browser and start the Java application.
- ▷ Wait until the Programme Manager is indicated.
- ▷ Define the bookmark.

Operation via VNC software

For operation of the furnace via VNC, a suitable VNC software must be installed on the computer.

- ▷ Start a VNC software (e. g. RealVNC or UltraVNC) on a computer in the network. The start window is displayed.
- ▷ In the *Server* entry line, enter the noted down IP address (for example "192.168.115.200").
- ▷ Click *OK*.
The login window is displayed.
- ▷ Enter the password in the *Password* entry line. The default password is "DEKEMA".
- ▷ Click *OK*.
The Programme Manager is displayed in the VNC window.
- ▷ Installing a fixed IP address for the furnace in the Programme Manager (see Chapter 5.4.3.).

5.4.3 Installing a fixed IP address for the furnace

Depending on the network settings, a different IP address is used for each restart of the unit. To avoid repeated finding and transferring of the IP address, a fixed IP address should be assigned to the furnace in the network.

- ▷ Have the network administrator assign an IP address, a subnet mask and a gateway address for the furnace.
- ▷ In the *Setup* menu, call up the *Netzwerk (Network)* function.

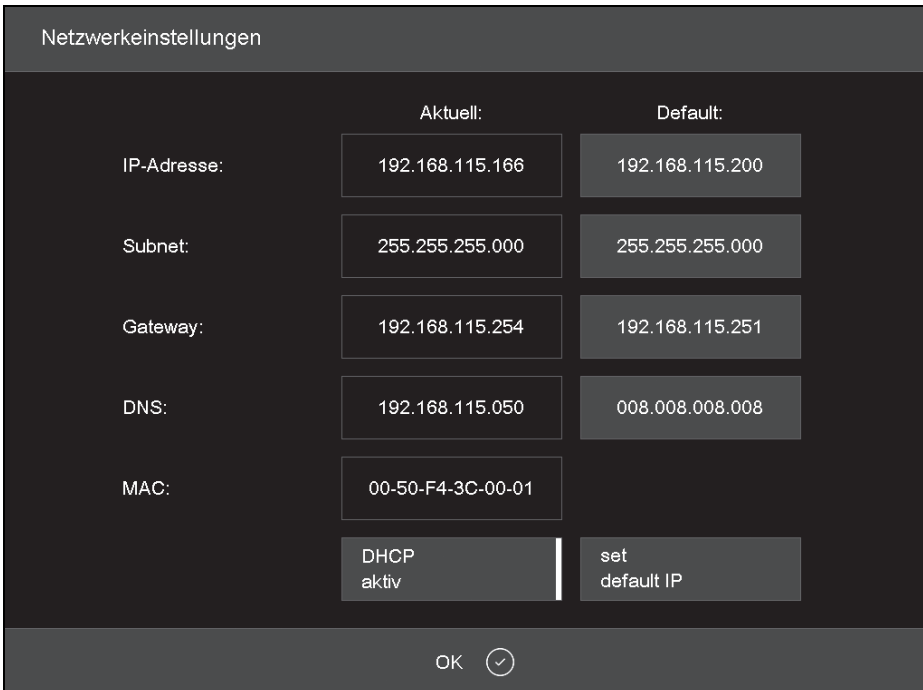
Fig. 6 Menu Setup > Network

- ▷ In *Default IP-Adresse (Default IP address)*, enter the IP address assigned to the furnace.
- ▷ In *Default Subnet Maske (Default Subnet Mask)*, enter the subnet mask assigned to the furnace.
- ▷ In *Gateway*, enter the gateway address assigned to the furnace.
- ▷ Tap on the *DHCP aktiv (DHCP active)* button to deactivate the DHCP function.

- ▷ Click *OK*.
The furnace now has a fixed IP and can be responded to via a browser or the VNC software under this address.

 The new IP address will not be activated until after restarting the furnace.

5.5 Connecting the furnace with an individual computer (optional)

- ▷ Make sure that the furnace is connected with the computer via a network cable.

 The factory-assigned IP address of the furnace is 192.168.115.200. If this IP address is no longer valid:

- ▷ Determine the current IP address of the furnace (see Chapter 5.4.1).

 If necessary, the IP address must be released in the computer's firewall (see operating system's manual).

- ▷ Adjust computer to correct IP address range as required.

5.5.1 Setting up the furnace operation

Operation via browser

 For operation of the furnace via Web interface, your browser must be Java-capable and Java must be activated.

- ▷ Start a Java-capable browser (e. g. Internet Explorer or Mozilla Firefox) on the computer.
- ▷ Enter the specified or determined IP address in the browser's address line (in the following example "http://192.168.115.200") and press *Enter*.

The Login window is displayed.

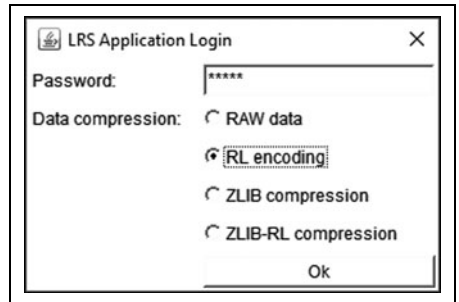
Fig. 7 Login window

- ▷ Enter the password in the *Password* entry line.
- ▷ Select the desired data compression for the data transfer between the furnace and the computer.

We recommend the standard setting *RL encoding*.

- ▷ Click *OK*.

The Programme Manager is displayed in the browser.

- ▷ In the *Setup* menu, call up the *Login Daten (Login Data)* function.

- ▷ Activate the functions *Start/Stop Funktionalität (Start/Stop Functionality)* and *Funktionstasten über Webserver (Function keys via Web server)*.

For comfortable calling up of the furnace control system, you can define a bookmark in the browser:

- ▷ Enter the IP address in the browser and start the Java application.
- ▷ Wait until the Programme Manager is indicated and define the bookmark.

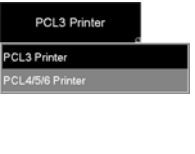
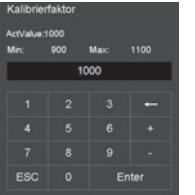
Operation via VNC software

For operation of the furnace via VNC, a suitable VNC software must be installed on the computer.

- ▷ Start a VNC software (e. g. RealVNC or UltraVNC) on the computer.
The start window is displayed.
- ▷ In the *Server* entry line, enter the IP address (in the following example “192.168.115.200”).
- ▷ Click *OK*.
The login window is displayed.
- ▷ Enter the password in the *Password* entry line.
The default password is “DEKEMA”.
- ▷ Click *OK*.
The Programme Manager is now displayed in the VNC window. See Chapter 5.3.

6 Operating elements

Name	Example	Description
Button		<p>A button is always associated with a function.</p> <ul style="list-style-type: none"> ▷ Tap on the button to execute the function.
On/Off button		<p>An On/Off button button is associated with a function. A bar and the colour indicate the status of the function:</p> <ul style="list-style-type: none"> _ Coloured: The function is activated. _ White: The function is deactivated. <ul style="list-style-type: none"> ▷ Tap on the button to change the status of the function.
Drop-down list		<p>Select the desired list entry as follows:</p> <ul style="list-style-type: none"> ▷ Tap on the entry. The complete selection list is opened. ▷ Tap on the desired entry. The list is closed and the selected list entry is indicated.
Entry field		<p>Enter/change the desired entry as follows:</p> <ul style="list-style-type: none"> ▷ Tap on the entry. The (alpha-)numerical keyboard window is displayed. The entry can be changed in it.
Numerical keyboard window		<p>The numerical screen keyboard is used for entering number values, e.g., firing parameters:</p> <ul style="list-style-type: none"> _ <i>ActValue</i> indicates the currently set value _ <i>Min</i> indicates the lowest allowed value _ <i>Max</i> indicates the highest allowed value
Alphanumerical keyboard window		<p>The alphanumerical screen keyboard is used for entering text, e.g., names of firing programmes:</p> <ul style="list-style-type: none"> ▷ Enter the desired text. ▷ Tap on <i>Enter</i> to finalise the entry ▷ Tap on <i>ESC</i> to cancel the entry

Tab. 2

Name	Example	Description
Keyboard window for classic editor		<p>▷ To edit the command line: Tap on Edit. The text is displayed within a blue frame.</p> <p>▷ Move the cursor to the desired location using the cursor buttons</p> <p>-or-</p> <p>▷ tap on the desired position in the graphical representation of the programme sequence.</p> <p>▷ Edit the command as desired.</p> <p>▷ To finalise the entry: Tap on <i>Enter</i>. The classic keyboard window is closed.</p>
	Swipe pad, lift	<p>To move the lift upward:</p> <p>▷ Swipe the swipe pad from the bottom upward.</p> <hr/> <p>To move the lift downward:</p> <p>▷ Swipe the swipe pad from the top downward.</p> <p>-or-</p> <p>▷ Press the <i>F4</i> key of the computer.</p> <hr/> <p>To stop the lift:</p> <p>▷ Press the swipe pad.</p> <hr/> <p>To move the lift to a defined position (top, closed/ centre/ bottom, opened):</p> <p>▷ Press the swipe pad at the top/in the centre/at the bottom.</p>
	On/Off button	<p>To end all current commands:</p> <p>▷ Briefly press the On/Off button.</p> <p>-or-</p> <p>▷ Press the <i>F9</i> key of the computer.</p> <hr/> <p>To switch to/out of standby:</p> <p>▷ Press the On/Off button for approx. 5 seconds.</p>

Tab. 2

7 Screen displays

The screen display for operation of the furnace is subdivided into different areas.

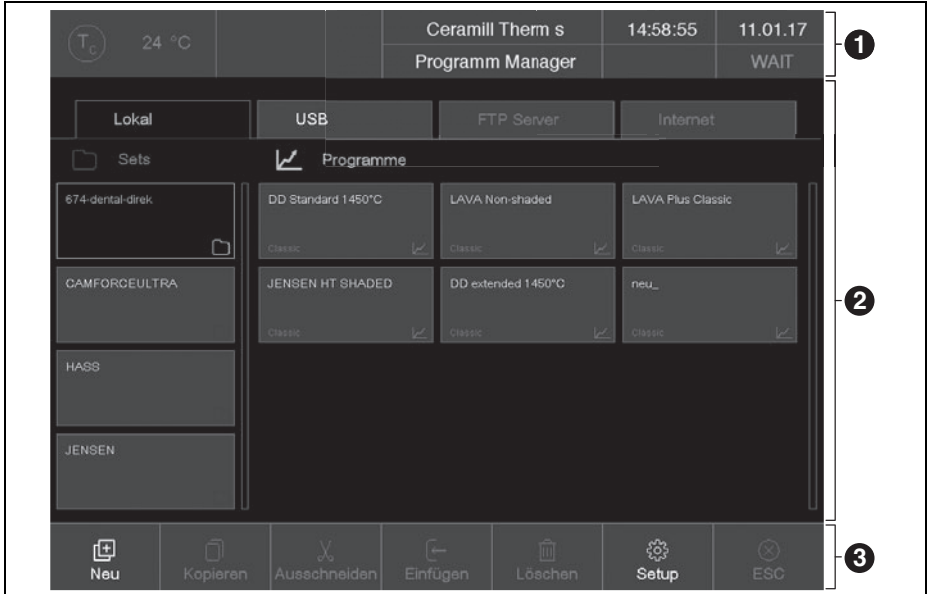
Fig. 8 Screen display (example: Programme Manager)

The information bar [1] displays all important parameters:

- _ The temperature in the firing chamber: *T_c*
- _ The furnace type or the firing number #. By tapping the button, you can switch between the displays.
- _ Time of the furnace control system
- _ The date of the furnace control system or - should a programme be selected - the presumable remaining time respectively the presumable end of the programme. By tapping the button, you can switch between the displays.
- _ The current menu (here: *Programm Manager (Programme Manager)*). A lock indicates a programme protected by a general code

_ The status of the furnace:

- *Wait*: The furnace is in waiting status
- *Run*: A programme is being executed (programme sequence)
- *Edit*: A programme is being edited (programme editor)

Display area [2] indicates the current menu of the furnace control system, e.g., the Programme Manager.

The toolbar [3] contains the buttons for controlling and operating the furnace. The names of the buttons depend on the selected menu.

7.1 Programme Manager

The Programme Manager is used for managing firing programmes and sets (such as copy, delete) and for programme navigation.

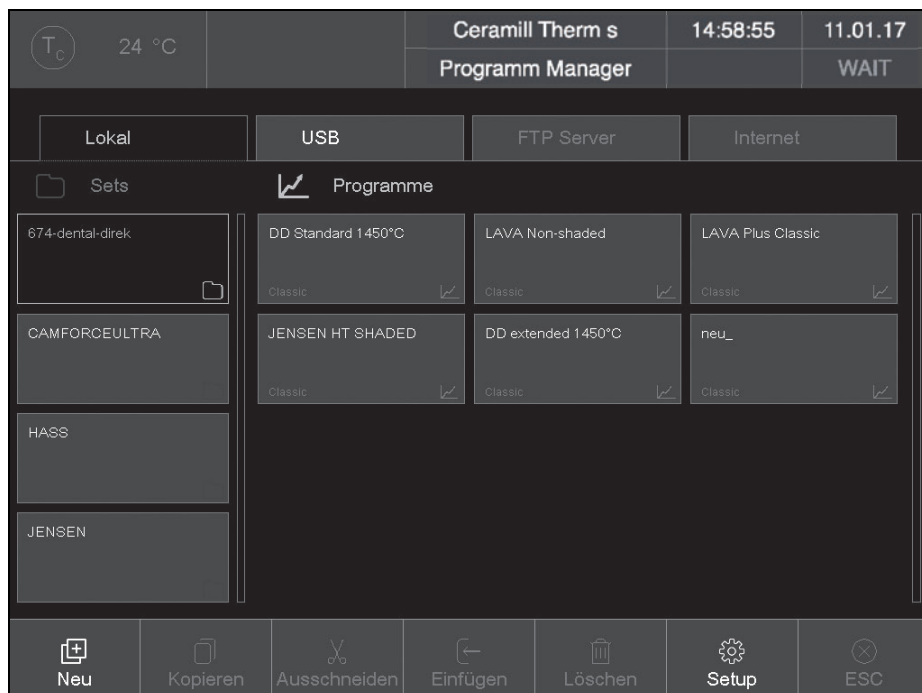
Fig. 9 Programme Manager (example)

7.1.1 Display range

The display range of the Programme Manager comprises the following elements:

- _ Storage area of programmes and sets:
 - *Lokal (Local)* displays the programmes and sets that are stored on the furnace
 - When a USB stick is inserted in the furnace interface, *USB* indicates the programmes and sets stored on this USB stick.
 - When the furnace is connected with the Internet, accessing the programmes and set on the DEKEMA server is possible *via Internet*.
 - When the furnace is connected with a FTP server, *FTP Server* indicates the programmes and sets stored on this FTP server.
- _ Set and programme list for the selected storage area:
 - The set list *Sets* displays all sets of the selected storage location. Sets are marked by a folder symbol  in the button.
 - The programme list *Programme (Programmes)* shows all firing programmes of the loaded set. Firing programmes are marked by a programme symbol  in the button. The programme list can contain folders for better clarity.

The scrollbar enables you to scroll through the respective list (swipe function).

Selecting programmes

To select a program:

- ▷ Tap on the desired storage area.
- ▷ Scroll in the set list until the desired set is displayed.
- ▷ Tap on the set to select it.
 - The respective button is displayed coloured and the corresponding firing programmes are indicated in the programme list.
- ▷ Scroll in the programme list until the desired firing programme is displayed.
- ▷ Tap on the firing programme to select it.
 - The programme indication is displayed.

Marking programmes

To mark a programme:

- ▷ Tap on a firing programme for approx. two seconds.
 - The programme is marked but not opened; the respective button is displayed coloured. Multiple selection of programmes is possible.

7.1.2 Toolbar

The buttons in the toolbar refer to the selected element *Sets* or *Programme (Programmes)*.

Programmes that are protected by the general code cannot be deleted or overwritten.

The **Neu (New)** button

To create a new set:

- ▷ Tap on *Sets*.
- ▷ Tap on *Neu (New)*.

The *Speichern unter (Save as)* window is displayed.

- ▷ Enter the name for the new set and tap on *Speichern (Save)*.

To create a new firing programme:

- ▷ Refer to Chapter 8.4.

To create a new folder:

- ▷ Tap on *Programme (Programmes)*.
- ▷ Tap on *Neu (New)*.

The *Programmauswahl (Programme selection)* window is displayed.

- ▷ Select *Ordner (Folder)*.
- ▷ Tap on the folder name.

The alphanumerical keyboard window opens.

- ▷ Enter the name for the new folder (40 characters max.) and tap on *OK*.

The new folder is saved/stored.

The **Kopieren (Copy)** button

To copy an element:

- ▷ Mark the desired element.
- ▷ Tap on *Kopieren (Copy)*.
- ▷ Navigate to the target destination and tap on the area where the copied element is to be inserted.

A blue bar shows the insertion point.

- ▷ Tap on *Insert*.

The **Ausschneiden (Cut)** button

To cut an element:

- ▷ Mark the desired element.
- ▷ Tap on *Ausschneiden (Cut)*.
- ▷ Navigate to the target destination and tap on the area where the cut element is to be inserted.

A blue bar shows the insertion point.

- ▷ Tap on *Einfügen (Insert)*.

The **Einfügen (Insert)** button

With the *Einfügen (Insert)* button, you can insert a cut or copied element at the insertion point.

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To delete an element:

- ▷ Mark the desired element.
- ▷ Tap on *Löschen (Delete)*.

A security prompt is displayed.

- ▷ When the element is to be deleted: Tap on *OK*.

The **Setup** button

The setup button opens the *Setup* menu (see Chapter 9).

The **ESC** button

In the Programme Manager, the *ESC* button has no function.

7.2 Programme indication

When a programme is selected in the Programme Manager, the programme indication with the firing parameters is displayed.

7.2.1 Classic programmes

The figure below shows a typical programme indication for classic programmes.

Fig. 10 Programme indication for classic programmes (example)

- 1 Information line
- 2 Graphical indication of the programme sequence
- 3 Command line with firing parameters
- 4 Toolbar

Information line [1]

The information line displays the programme name.

The two sequence buttons < and > left and right of the information line scroll to the next respectively the previous programme. Alternatively, you can also swipe left or right across the screen for this.

Graphical indication [2]

The graphical indication shows the programme sequence and contains the following elements:

- _ **Vertical axis:** Indicates the temperature.
- _ **Horizontal axis:** Indicates the time.
- _ **Cursor:** Indicates the current position (time) in the programme sequence.
- _ **Diagram:** Displays the programme sequence as a coloured curve. The chronological sequence of the temperature with respect to the firing object is represented. Time intervals in the graphic are backed in colour.
- _ **<, >:** The two sequence buttons aside of the graphic move the cursor through the individual programme steps.

 Double-tapping in the graphical representation opens the programme editor in order to change the firing parameters.

Command line [3]

The command line displays the individual programme steps with the firing parameters in the “CVTL” syntax. A tap on the command line opens the programme editor for changing firing parameters.

 Tapping in the command line opens the programme editor in order to actually change the firing parameters.

Toolbar [4]

The toolbar contains various buttons:

- _ *Speichern (Save)* saves/stores the programme.
- _ *Edit* opens the programme editor in order to change the firing parameters.
- _ *GO TO* opens the *GO TO Auswahl (GO TO selection)* window. In this window, a different programme can be selected within the set.

Fig. 11 GO TO Auswahl (GO TO selection) window

- _ *Bilder (Images)* opens a window for selection of an image.
- _ *Start* starts the programme. After starting, the *Start* button is automatically replaced through the *Stop* button.
- _ *Stop* prematurely ends a running programme. The heating is switched off and the furnace cools down. The firing run is not counted, meaning, the firing number does not increase. After stopping, the *Stop* button is automatically replaced through the *Start* button.
- _ *ESC* ends the programme indication and opens the Programme Manager.

7.2.2 Images

Files can be stored in the memory of the furnace electronics:

- _ Pictures, for example, in order to compare firing results with photos from sample firings
- _ PDF files, as an example, for readily access to this User Manual.

The screen page offers simplified functionality for selecting and opening files, as known for example from Windows Explorer. To import or transfer the files, you can use the FTP functionality or a USB stick.

7.3 Programme sequence

A firing programme is started from the programme indication by pressing *Start*. The programme steps are executed according to the specified sequence in the entry mask..

A running programme can be prematurely ended with the *Stop* button. In this case, all commands are terminated and the heating is switched off. The furnace cools down. The firing run is not counted, meaning, the firing number does not increase.

After the end of the programme, the QM information is displayed (if enabled).

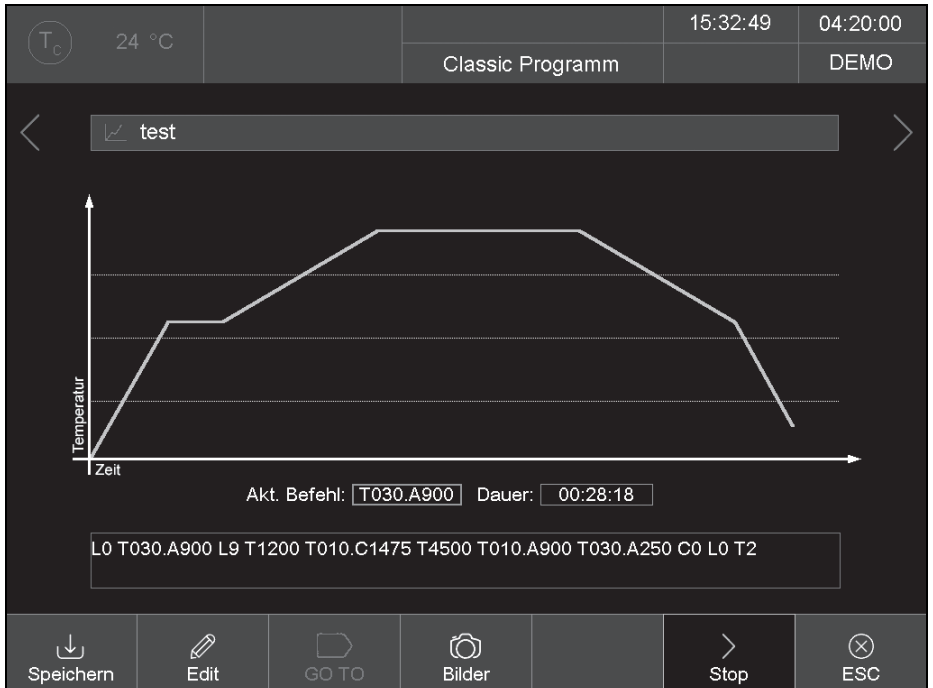
Fig. 12 Programme indication for classic programmes (example)

The example in Fig. 12 leads to the following firing procedure:

1. **LO**: Programme start with firing chamber open, meaning, the lift moves down if applicable.

By default, the firing procedure will not be executed until after the correct temperature for starting the selected programme has been reached. Otherwise, the Tc indicator will flash and a message will be displayed that the furnace is too hot or still heating up.

1. **T030.A900**: Adjusts the object temperature (900 °C) with a ramp of 30 °C/min.
The lift partially moves up, while the firing chamber is being heated up at the same time.
2. **L9**: Close (the lift moves upward).
3. **T1200**: Waiting period 1200 seconds.
4. **T010.C1475**: Adjusts the firing temperature (1475 °C) with a ramp of 10 °C/min.
5. **T4500**: Waiting period 4500 seconds.
6. **T010.A900**: Adjusts the object temperature (900 °C) with a ramp of 10 °C/min.
The lift partially moves down, while the firing chamber cools off at the same time.
7. **T030.A250**: Adjusts the object temperature (250 °C) with a ramp of 10 °C/min.
The lift partially moves down, while the firing chamber cools off at the same time.
8. **CO**: Switches off the heating.
9. **LO**: The lift moves completely down.
10. **T2**: The acoustic signal indicating the end of the programme sounds. The firing number is automatically increased by one.

8 Programme editor

8.1 CVTL commands

The firing programmes are programmed using CVTL commands. These commands have the following meaning:

Command	Meaning
A	defines the temperature in °C on the insulation table (A utodry)
C	defines the temperature in °C in the firing chamber (degree C elsius)
L	L defines the lift position (L ift): _ L0 = Lift completely down _ L9 = Lift completely up
T and an at least two digit number	defines the time in seconds (T ime)
T and a single digit number	defines the number of signal tones (beeps)
TOxx.Cxxxx	defines a temperature ramp via temperature gradient: _ T defines the ramp in °C/min, _ C defines the final temperature in °C in the firing chamber
T(x)xx.Cxxxx, (x) <> 0	defines a temperature ramp: _ T defines the time in seconds, _ C defines the final temperature in °C in the firing chamber
TOxx.Axxxx	defines a temperature ramp via temperature gradient: _ T defines the ramp in °C/min _ A defines the final temperature in °C on the insulation table
T(x)xx.Axxxx, (x) <> 0	defines a temperature ramp: _ T defines the time in seconds _ A defines the final temperature in °C on the insulation table
T(x)xx.Lx	defines a lift ramp: _ T defines the time in seconds _ L defines the final lift position

Tab. 3 CVTL commands

For better clarity, two successive commands can be separated with a space.

On confirmation of the entry, the syntax is checked. If the entered command line contains an error, a warning message is displayed. The command line must be corrected.

8.2 Calling up the programme editor

- ▷ Select the desired firing programme in the Programme Manager.

The programme indication is displayed.

- ▷ Tap on *Edit*.

A coloured button indicates that the programme editor is active. Additionally, the *Edit* status is displayed in the information line.

- ▷ Tap on the command line below the programme indication.

The keyboard window for classic editor is displayed

The entry of a new command line is expected by default.

Commands are overwritten by default. To enter new commands:

- ▷ Tap on *Insert*.

The command line shows the command sequence in the order of its execution. Programming with the CVTL commands being used requires practice and understanding of the ceramics and firing techniques.

8.3 Cancelling the programme editor

- ▷ Tap on *Speichern (Save)* save/store the changes.

- ▷ Tap on *Edit*.

The coloured button goes out and indicates that the programme editor is inactive. The *Edit* status is no longer displayed in the information line.

8.4 Creating a new programme

To create a new programme:

- ▷ Open the Programme Manager.
- ▷ Tap on the desired storage area.

- ▷ Tap on the set where the new programme is to be saved/stored.

- ▷ Tap on *Programme (Programmes)*.

The programme list is active now.

- ▷ Tap on *Neu (New)*.

The *Programmauswahl (Programme selection)* window opens.

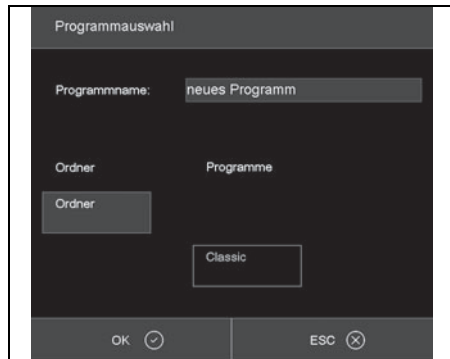
Fig. 13 Programme selection window

- ▷ Tap on *Classic*.

- ▷ Tap on the programme name.

The alphanumeric keyboard window opens.

- ▷ Enter a programme name.

- ▷ Tap on *OK*.

The programme editor is opened.

- ▷ Edit the new programme.

- ▷ Tap on *Speichern (Save)*.

- ▷ Tap on *ESC*.

The programme editor is closed and the Programme Manager is displayed.

The programme name can be changed later at any time:

- ▷ In the Programme Editor, tap into the field for the programme name.

The alphanumeric keyboard window opens.

9 Setup menu

The following settings can be made in the *Setup* menu:

- _ System control
- _ Global furnace settings
- _ Maintenance programmes
- _ Data recording

The *Setup* menu is called up by tapping on the *Setup* button in the Programme Manager.

 The menu may be customer-specifically blocked by the supplier.

Fig. 14 Setup menu

Toolbar buttons:

- _ *Service* calls up the *Service* menu.
The *Service* menu is only available with level 2 access right.
- _ *ESC* takes you back to the Programme Manager.

The table below lists the available function in the *Setup* menu along with a brief description.

Function	Description
Einstellungen (Settings)	
<i>Einheiten (Units)</i>	For selection of the temperature units.
<i>Sprache (Language)</i>	For selection of the user interface language.
<i>Datum/Uhrzeit (Date/Time)</i>	For setting the date and time in the furnace control system.
<i>Bildschirmeinstellung (Screen setting)</i>	For setting of the time, after which the screensaver is activated without an entry.
<i>Netzwerk (Network)</i>	For adjustment of the furnace network settings.
<i>Generalcode (General code)</i>	For adjustment of the furnace access levels.
<i>Ofenkennung (Furnace ID)</i>	For entry of the furnace ID.
<i>Akustisches Signal (Acoustic signal)</i>	For setting the acoustic signal at programme end.
<i>Heizungseinstellungen (Heating settings)</i>	For definition of the maximum values for temperature rise rate and power consumption.
<i>Lifteinstellungen (Lift settings)</i>	Determines whether a hazard warning is displayed before the lift opens while the firing chamber is hot.
Programme (Programmes)	
<i>Trocknen (Drying)</i>	Performs a programme run for drying the furnace.
<i>Checkprogramm (Check programme)</i>	Performs a programme to check the operability of the most important assemblies.
<i>Oxidation</i>	Performs the oxidation programme.
<i>Temperaturkalibrierung (Temperature calibration)</i>	Performs a temperature calibration and adjusts the CCF value (calibrations factor).
Kommunikation (Communication)	
<i>Drucker (Printer)</i>	Specifies a printer for the furnace.
<i>Qualitätsmanagement (Quality management)</i>	Determines whether a QM protocol is displayed after the end of a programme sequence.
<i>Diagnosedaten (Diagnosis data)</i>	Displays the diagnostic data that were created from the check programme.
<i>Login Daten (Login data)</i>	Defines the access rights to the furnace via the web interface.
<i>Backup / Recovery</i>	Backup and recovery of the data of the furnace control system.
<i>FTP Server Login</i>	Defines the furnace access to a FTP server.

Tab. 4

9.1 Settings

9.1.1 Network

▷ In the *Setup* menu, tap on the *Netzwerk (Network)* button.

The window with the network settings of the furnace is displayed. The current network addresses are displayed.

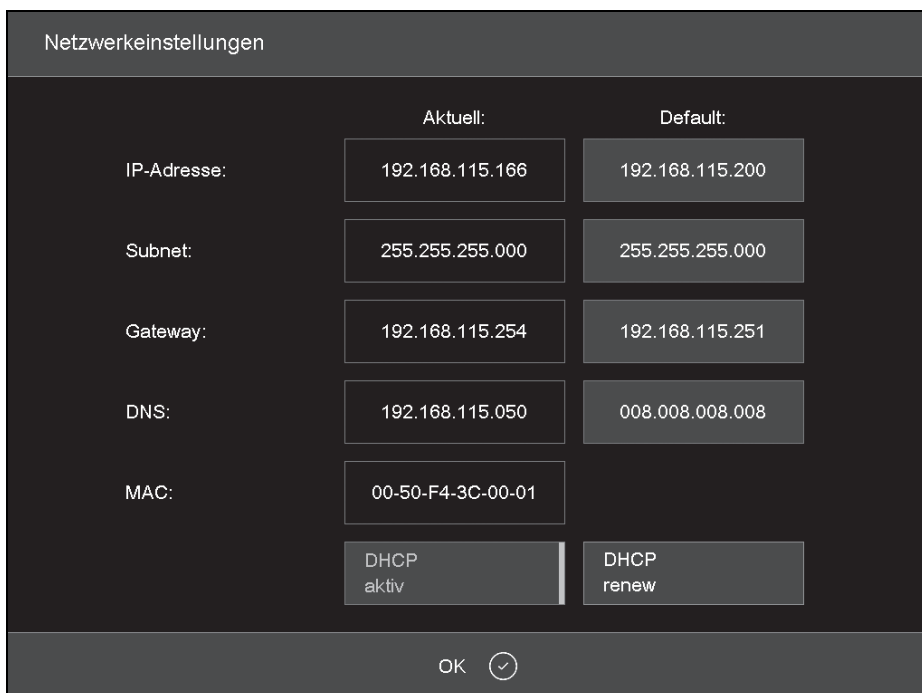
Fig. 15 Setup > Network

Buttons:

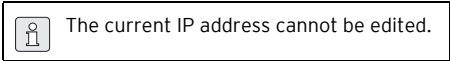
- *OK* accepts the settings and returns to the *Setup* menu.
The changed settings will not be activated until after restarting the furnace.
- *DHCP aktiv (DHCP active)* allows the furnace to have an IP address assigned by your network

server. A coloured button indicates that the DHCP function is active.

- *DHCP renew* (only when DHCP is active) allows the furnace to have a **new** IP address assigned automatically by your network server.
- *set default IP* (only when DHCP is not active) assigns the default IP address to the furnace.

To change an IP address:

- ▷ Tap on the respective address.
The alphanumerical keyboard window is displayed.
- ▷ In the keyboard window, tap on *Edit*.
- ▷ Enter the address and confirm with *Enter*.

9.1.2 General code

This function allows you to protect (lock) firing programmes against being changed.

- _ With the general code, only the programmes on the *Local* memory of the furnace control system can be protected.
- _ With the general code, all programmes in a set are protected or released. A different marking of individual programmes is not possible.
- _ Protected programmes can not be edited and changed. The programme editor cannot be called up for these programmes.
- _ Protected programmes cannot be deleted (Programme Manager).
- _ The functionality in the *Setup* menu is heavily limited.

Sets protected by the general code are marked with a lock in the button in the Programme Manager. In the programme indication, a lock in the information bar indicates a protected programme.

There are two different levels of protection:

- _ Protected programmes
Protected programmes cannot be changed or deleted. They are marked by the “” symbol.
- _ Fixed programmes
The firing parameters of fixed programmes can be adapted for a one-time execution, but not saved/stored. After executing such a programme, the Programme Manager is automatically displayed. Fixed programmes are marked “ Fix”.

To activate the general code:

- ▷ In the *Setup* menu, tap on the *General code* button.
The alphanumerical keyboard window is displayed.
- ▷ Enter a password for the general code.
The *General code* window is displayed.

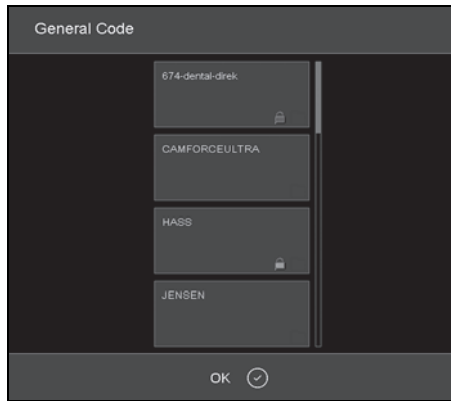
Fig. 16 General code window (example)

To protect the programmes in a set:

- ▷ Tap on the respective set.
The symbol in the button indicates that all programmes in this set are protected.

- ▷ Tap once more on this set.
The word *Fix* is displayed next to the  symbol and indicates that all programmes in this set are treated as fixed programmes.
- ▷ Tap once more on this set.
The  symbol disappears. The programmes in this set are no longer protected.
- ▷ Protect all desired sets in this manner.
- ▷ Tap on *OK*.
The *Setup* menu is displayed.
- ▷ Sets protected by the general code are marked with a lock in the button in the Programme Manager. In the programme indication, a lock in the information bar indicates a protected programme.
- ▷ To activate the general code:
- ▷ In the *Setup* menu, tap on the *General code* button.
The alphanumerical keyboard window is displayed.
- ▷ Enter the password for the general code.
The general code is deactivated.

To change protected sets:

- ▷ Deactivate the general code.
- ▷ Change programmes or add sets.
- ▷ Activate the general code and protect the sets as desired.

9.2 Programmes

9.2.1 Drying

The drying programme is used for drying the firing chamber, for example when the furnace was subject to high humidity for a prolonged period.

- ▷ In the *Setup* menu, tap on the *Trocknen (Drying)* button.
The window for starting the drying programme is displayed.

Fig. 17 Setup > Drying

 *ESC* closes the window and switches back to the *Setup* menu.

- ▷ Mount the firing table.
- ▷ Tap on *Start*.
The *Trocknen (Drying)* window is displayed as long as the programme is running (approx. 90 minutes).

The programme run can be cancelled at any time:

- ▷ Tap on the *Stop* button.
A security prompt - if you really want to cancel the programme - is displayed:
 - *OK* confirms the prompt and cancels the programme.
 - *ESC* withdraws the decision and continues the programme.

9.2.2 Check programme

▷ See Chapter 11.2.2 on Page 83.

9.2.3 Temperature calibration

▷ See Chapter 11.2.1 on Page 81.

9.3 Communication

9.3.1 Quality management

▷ In the *Setup* menu, tap on the *Qualitäts-Management (Quality management)* button.

The quality management window is displayed.

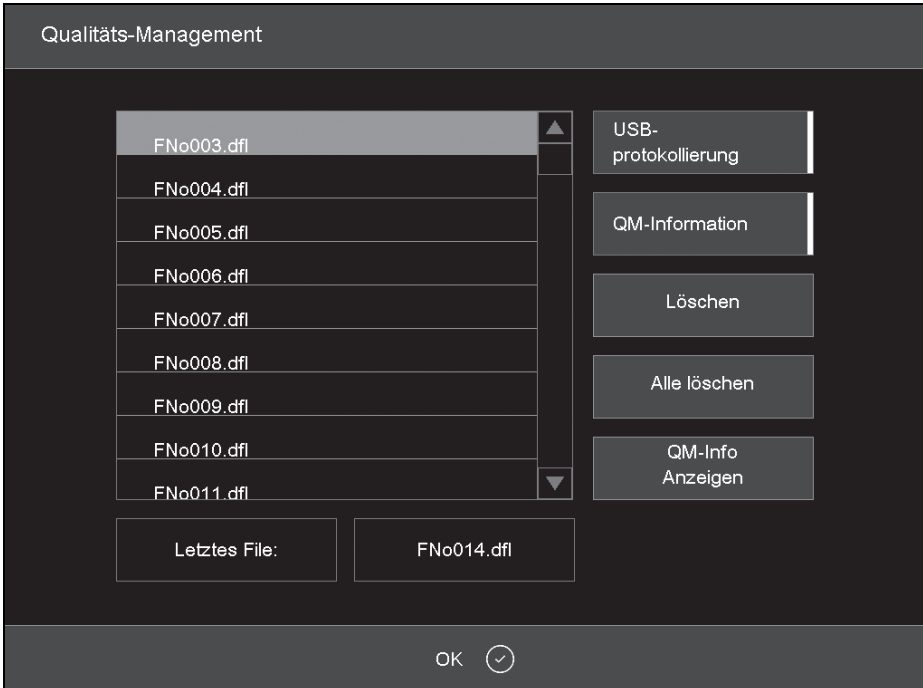
Fig. 18 Setup > Quality management

The furnace logs the last ten completed firing programmes. For more than ten completed firing programmes, these QM files are overwritten on a rolling basis.

_ *OK* accepts the settings and returns to the *Setup* menu.

_ *Letztes file (Last file)* displays the QM file name of the last completed firing programme.

_ *USB-Protokollierung (USB logging)* activates/deactivates the data recording on an external USB stick. A coloured button indicates that this function is active. Through this, the QM files are additionally updated on this USB stick (USB logging).

- _ *QM information* activates/deactivates the display of the QM information after each completed firing programme. A coloured button indicates that this function is active.
 - Not active: QM information is not displayed, but logged.
 - Active: After a completed firing programme, the corresponding detailed firing information are indicated on the display.
- _ *Löschen (Delete)* deletes the selected QM file.
- _ *Alle löschen (Delete all)* deletes all QM files.
- _ *QM-Info anzeigen (Display QM info)* displays the selected QM file.

9.3.2 Backup / Recovery

 For reasons of data security, we recommend that you back up all data regularly.

- ▷ In the *Setup* menu, tap on *Backup/Recovery*.
- ▷ The window for backing up and recovering data is displayed.

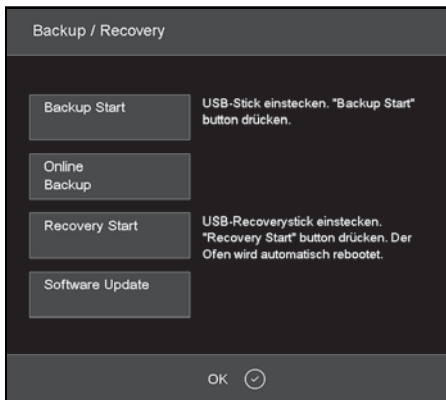
Fig. 19 Setup > Backup / Recovery

Buttons:

- _ *Backup Start* starts backing up the data on a USB stick.

- _ *Online Backup* starts backing up the data via the network connection (online).
- _ *Recovery Start* starts recovering the data on a USB stick.
- _ *Software Update* starts an online update of the furnace software.
- _ *ESC* closes the window and switches back to the *Setup* menu.

To back up data on a USB stick:

- ▷ Insert a USB stick with sufficient memory into a USB port on the back side of the furnace.

 Any existing data on the USB stick will be deleted at the beginning of the backup.

- ▷ Tap on *Backup Start*.
All furnace data are written on the USB stick.
- ▷ Wait until an acknowledgement tone sounds.
It signals the completion of the backup process.
- ▷ Remove the USB stick.
The furnace data are backed up.

To recover data from a USB stick or perform a software update:

- ▷ Insert a USB stick with data backup/software update into a USB port on the back side of the furnace.
- ▷ Tap on *Recovery Start*.
All furnace data will be restored and/or the software is updated.
- ▷ Wait until an acknowledgement tone sounds.
It signals the completion of the process.
- ▷ Remove the USB stick.
All furnace data are restored and/or the software is updated.

To back up data via online backup:

 An online backup requires an active Internet connection.
Only backup and software update are possible online; not the recovery of data.

- ▷ Tap on *Online Backup*.
All furnace data (operating system, software, set parameters and firing programmes) are stored on the server.
- ▷ Wait until the backup process is completed. This can take several minutes.
The furnace data are backed up.

To update the software online:

 An online update of the furnace software requires an active Internet connection.

- ▷ Tap on *Software Update*.
A message indicates that the furnace is looking for a new version of the software. When a new version of the furnace software is available, a prompt to download the new version is displayed.
- ▷ Tap on *OK*.
- ▷ Wait until the update process is complete and the furnace reboots. This can take several minutes.
The software is updated.

9.3.3 FTP server login

The FTP functionality provides access to an externally connected FTP server. As an example, this allows for firing programmes, image folders and files to be linked to the furnace.

 If necessary, request the login data for the FTP server from the IP administrator.

- ▷ In the *Setup* menu, tap on the *FTP Server Login* button.
The login window for a FTP server is displayed.

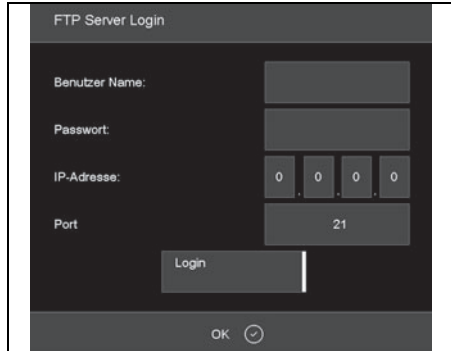
Fig. 20 Setup > FTP server login

- Buttons:
- _ *ESC* closes the window and switches back to the *Setup* menu.

- To establish a connection to the FTP server:
- ▷ Tap on the *Benutzer Name (User name)*, *Passwort (Password)*, *IP-Adresse (IP address)* and *Port* entry fields, and enter the respective values.
 - ▷ Tap on *Login*.
The furnace is connected with the respective FTP server. A coloured *Login* button indicates that the connection is active.

- To end the connection:
- ▷ Tap on *Login*.
The furnace is disconnected from the respective FTP server. The coloured *Login* button goes out.

10 Operation

10.1 Switching on and off

10.1.1 Preparations

- ▷ Check if the mains cable of the sintering furnace is plugged into a socket outlet with earthing contact.
- ▷ Check if the furnace is connected via a web interface with a single computer or with a network (optional).

10.1.2 Switching on

Cold start:

- ▷ Switch the furnace on via the mains switch on the back side.

The furnace control system boots up (approx. one minute). Afterwards, an internal check of the furnace hardware and software takes place and the lift automatically moves to the bottom position.

Start from idle state:

- ▷ Switch the furnace on with the On/Off switch next to the touchscreen.
The furnace is immediately ready for operation.

- ▷ Check, if the insulation table is seated on the lift plate. If this is not the case:
 - Switch the furnace off.
 - Mount the insulation table.
 - Switch the furnace on again.

10.1.3 Switching off

- ▷ Shut the firing chamber.
- ▷ Switch the furnace into the idle state with the On/Off switch next to the touchscreen.

-or-

- ▷ Switch the furnace off via the mains switch on the back side.

10.2 Placing on and removing a firing object

10.2.1 Placing on a firing object

NOTE:

Avoid damaging the firing object and the heat-insulation material!

- ▷ When positioning the firing object, pay attention that no parts project over the upper brim of the sintering plate.

- ▷ Place on the wave-formed sintering ring
- ▷ Mount the sintering plate.
- ▷ Place the firing object on the sintering plate.
- ▷ If required: Place on a second sintering ring with sintering plate and the firing objects.

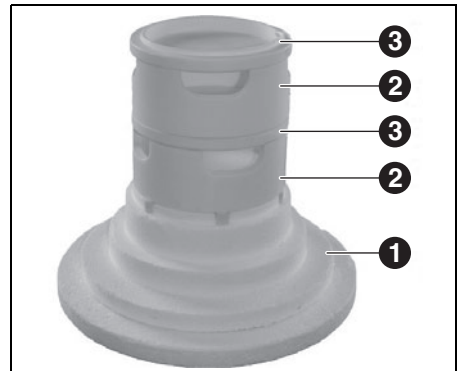
Fig. 21 Maximum stacking of the sintering rings

- 1 Insulation table
- 2 Sintering ring
- 3 Sinter plate

10.2.2 Removing the firing object

- ▷ Remove the firing object using suitable aids (tweezers, safety gloves).

10.3 Selecting and starting programmes

After switching on the furnace, the Programme Manager is automatically displayed on the touch-screen.

▷ Select the desired programme.

The firing procedure of the selected programme is displayed.

▷ In the programme indication, tap on *Start*.

The programme is automatically run.

NOTE:

During the run of a programme, the lift automatically moves up and down.

- ▷ Make sure that no objects obstruct the lift's motion.
- ▷ Do not attempt to manually stop or accelerate the lift.

During the programme run, the programme can be cancelled by tapping on *Stop*.

10.4 Programme end

For chamber temperatures $T > 999\text{ °C}$, the firing chamber generally remains shut. It cannot be opened until the chamber temperature is below 999 °C .

- _ After a regular programme ends, an acoustic signal sounds.
- _ The lift moves down.
- _ The firing number is automatically increased by one.

During the programme run, the programme can be cancelled by tapping on *Stop*. In this case, all commands are terminated and the heating is switched off. The furnace cools down. The firing run is not counted, meaning, the firing number does not increase.

10.5 After a power failure

After a power failure, the control system remembers the condition of the furnace at the time of interruption for up to two hours.

When the power supply is reactivated within these two hours, the furnace will restore the condition at the time of the interruption. If the temperature drop in the firing chamber does not exceed 150 °C , the firing programme being executed with continue.

Otherwise, the furnace will prompt the inquiry, if the programme is to be continued.

10.6 Access rights

Three levels are defined for accessing the furnaces' user interface. The access levels 1 and 2 are password-protected. The passwords are defined in the *Service* menu via the *Passwörter (Passwords)* function.

The access rights should be assigned to the operators of the furnace according to their knowledge, experience and responsibility.

By default, the furnace starts in the *Bedienung (Operation)* menu, with *General code* one level higher (level 0).

Access level	Description
General code (level 0)	Limited access to the <i>Setup</i> menu. Blocked programmes can be selected and started. Programmes not blocked can be edited and created. Resetting the general code
Operation (level 1)	All programmes can be selected and started. All programmes can be edited and created. Access to the <i>Setup</i> menu.
Service (level 2)	Access to the <i>Service</i> menu.

Tab. 5

11 Cleaning and Maintenance

 **WARNING:**
Possible damage when using non-original spare parts.
Malfunctions may occur when using non-original parts. Under certain circumstances, this can result in serious injury or substantial property damage.
▷ Use only original spare parts.

 **WARNING:**
Improper repairs.
▷ Have technical interventions carried out only through authorised specialised personnel.

11.1 Cleaning

 **DANGER:**
Risk of electric shock!
▷ Prior to any cleaning, pull the mains plug.

 **NOTE:**
Possible damage to the furnace!
▷ Do not use a high-pressure cleaner or compressed air for the cleaning!
▷ Use only commercially available, water-based or non-flammable, solvent-free cleaning agents for the cleaning.

Prior to cleaning:

- ▷ Allow the furnace to cool down to the room temperature.

- ▷ Clean housing and operating/control elements with a soft cloth and a mild cleaning agent. Do not use abrasive cleaners, scouring pads or solvents such as alcohol, petrol, or acetone.
- ▷ To ensure proper sealing, always keep the seal ring of the lift plate free of dust. For this, use a brush with a receptacle (ceramic dust) or a moist cloth.

- ▷ In the *Setup* menu, tap on the *Temperaturkalibrierung (Temperature calibration)* button. The temperature calibration window is displayed.

11.2 Maintenance

11.2.1 Temperature calibration

 For new units, we recommend performing an initial temperature calibration after approx. 6 months.

Required aids:

- _ Calibration set with PTC ring and calibration chart.
- _ Sliding calipers with a measuring accuracy exceeding $\pm 0.01 \text{ mm}$ ($= 10 \text{ }\mu\text{m}$), in order to ensure precise evaluation and temperature calibration.

 Prior to calibration, the furnace must be thoroughly heated through.

- ▷ Place the PTC ring centred on the insulation table.

Fig. 22 *Setup* > *Temperature calibration*

 *ESC* closes the window and switches back to the *Setup* menu.

- ▷ Tap on *Start*.
The programme for temperature calibration is started. It takes approx. 180 minutes.

The programme run can be cancelled at any time:

- ▷ Tap on the *Stop* button.
A security prompt - if you really want to cancel the programme - is displayed:
 - *OK* confirms the prompt and cancels the programme.
 - *ESC* withdraws the decision and continues the programme.

After the programme has run off:

- ▷ Allow the sintered PTC ring to cool down to room temperature on the insulation table.

- ▷ Remove the PTC ring from the insulation table using suitable aids (tweezers, safety gloves) and place it on a clean and level surface.
- ▷ Using the sliding calipers, measure the diameter D between the two opposing middle numbers of the three-digit batch number to an accuracy of ± 0.01 mm.

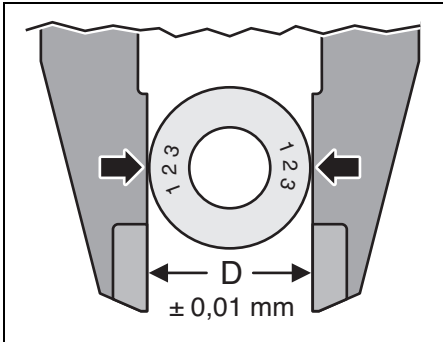
Fig. 23 Determining the diameter

- ▷ Compare the measured diameter D with the values of the ring-specific calibration chart. The respective value in the chart provides the actual temperature in the furnace.

If the temperature determined in this manner varies from the 1500 °C set-temperature, then the CCF value must be changed.

Determined

temperature	Corrective measure
> 1500 °C	Reduce CCF
= 1500 °C	No corrective measure required
< 1500 °C	Increase CCF

Tab. 6

To change the CCF value:

- ▷ In the *Setup* menu, tap on the *Temperaturkalibrierung* (*Temperature calibration*) button. The temperature calibration window is displayed. The set CCF value is displayed.

Fig. 24 Setup > Temperature calibration

- ▷ Tap on the CCF value. The numerical keyboard window is displayed.
- ▷ Enter the desired CCF value:
 - An increase of the CCF value by one point corresponds with an increase of the firing temperature by 1 °C.
 - An decrease of the CCF value by one point corresponds with an decrease of the firing temperature by 1 °C.
- ▷ Tap on *Enter* to save/store the new CCF value.

In case of large or sudden variations to the original CCF value:

- ▷ Check if there is a technical malfunction.
- ▷ Notify service partner.

11.2.2 Check programme

The furnace has an internal programme with which the operability of the most important components are checked. The checks include:

- _ The detection of aged components, such as wear of heating elements Correction of malfunctions.

After approx. 500 firing runs, the unit displays a recommendation to start the check programme.

▷ In the *Setup* menu, tap on *Check programme*.

The *Check programme* window is displayed.

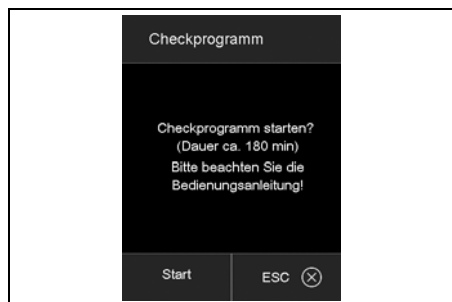
Fig. 25 *Setup* > *Check programme*

▷ Tap on *Start*.

The check programme is started. It takes about 180 minutes, depending on the set maximum increase rate.

On completion of the check programme, the results are saved/stored in a diagnosis file. If these diagnosis data show that the firing procedure of a programme presumably cannot be worked off correctly, then the furnace will automatically prompt this information when starting a firing programme. As an example, this will be the case when the heating element is already excessively worn and a high firing temperature needs to be reached quickly.

Error message	Corrective measure
KEINE FEHLER GEFUNDEN (NO ERROR FOUND)	No corrective measure required
HEIZKÖRPER WECHSELN (REPLACE HEATING ELEMENT)	Contact service partner

Tab. 7 *Error messages upon completion of the check programme*

12 Storage

- ▷ Switch furnace off and allow it to cool down to room temperature.

NOTE:

Several hours are required for the furnace to cool down.

- ▷ Switch the furnace on again.
- ▷ Move insulation table up in order to shut the firing chamber.
- ▷ Switch the furnace off.
- ▷ Pull mains plug of power supply from the mains socket outlet.

The furnace may only be stored in the original packaging. In this, the packaging symbols must be observed and the storage conditions must be met in accordance with Chapter 15.

13 Malfunctions, Repairs and Warranty

13.1 Malfunctions

- ▷ In case of malfunctions, contact the manufacturer's Technical Service.

13.2 Repairs

Repairs may be carried out only by trained, specialised personnel.

13.3 Warranty

The warranty complies with the statutory provisions. For more information, please refer to our General Terms of Business.

Each furnace is subject to a factory test run on behalf of the manufacturer. The test run may cause discolouration to the furnace compartment or sheet-metal parts.

14 Environmental Protection

Packaging

In terms of packaging, Amann Girrbach participates in country-specific recycling systems, which ensure optimal recycling.

All packaging materials used are environmentally friendly and recyclable.

Waste equipment

Waste equipment contain materials that should be sorted for environmental-friendly recycling.

- ▷ After the service life has elapsed, dispose of the machine via the public disposal systems in an environmentally-friendly manner.

15 Technical Data and Accessories/Spare Parts

 Subject to changes.

Technical data

	Unit	Value
Art. No.	–	178390
Dimensions (W × H × D)	mm	383 × 780 × 343
Usable size of the firing chamber (Ø × H)	mm	60 × 50
Weight	kg	23
Electrical connections	V/Hz	230/50-60
Max. output	W	2000
Max. temperature	°C	1560
Fuse (time-delayed)	A	2 × 16
Degree of protection	–	IP20
Thermal protection class according to DIN EN 60519-2	–	Class 0
Ambient conditions in operation:		
_ Temperature	°C	+15 ... +35
_ Humidity, non-condensing	%	60
_ Maximum operating elevation (standard elevation zero)	m	2000
Ambient conditions out of operation; storage:		
_ Temperature	°C	50
_ Humidity, non-condensing	%	80
Ambient conditions during transportation:		
_ Temperature	°C	+10 ... +55
_ Humidity, non-condensing	%	80
Interfaces	–	USB/Ethernet

Tab. 8

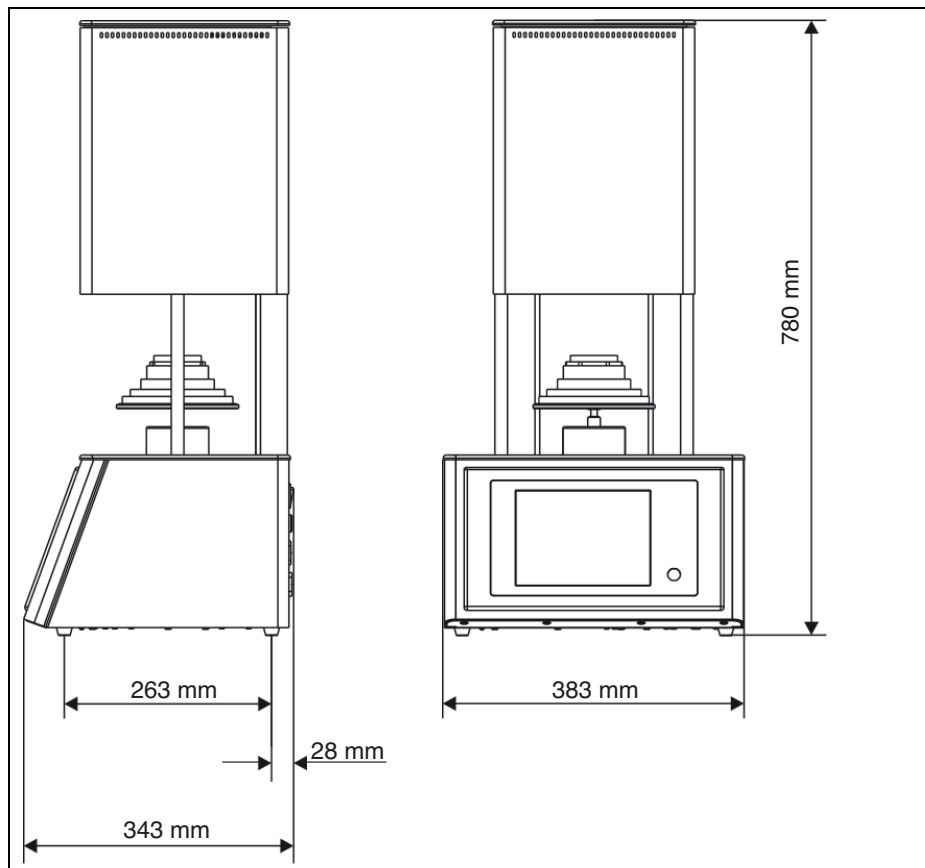
Fig. 26 Dimensions

Accessories/Spare parts

Art. No.	Designation
178391	Heating element with cable set 664/674
178392	Thermoelement 664/674 Vers.2 90°
178393	Firing table, ring d = 65 mm, 28 mm, 664

Art. No.	Designation
178394	Firing table, base, trough d = 65 mm, 664
178395	Insulation table, cpl. 664
178396	Fuses 5x20 16A 10 pce.

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Dental-Keramiköfen GmbH
Amann Girrbach GmbH